





RIGID BONDING SOLUTIONS FOR MEDICAL DEVICE ASSEMBLY



HIGH STRENGTH & DURABILITY

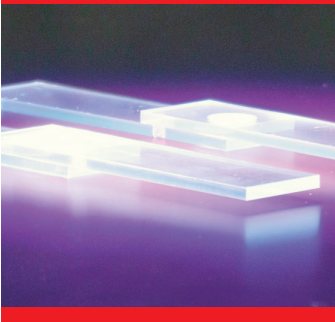


LOW VISCOSITY

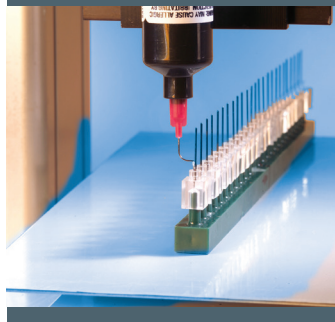

LED CURABLE (FAST & TACK-FREE)



HIGHLY FLUORESCENT

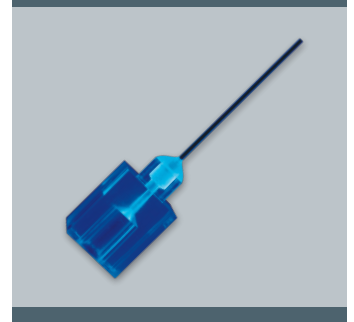

To learn more about our Tough and Flexible Adhesives, go to na.henkel-adhesives.com/medicaladhesives or call 1-800-LOCTITE




Henkel Adhesive Technologies

LOCTITE® AA 3961™ and AA 3963™ are the latest light cure innovations from LOCTITE. They have been tested to LOCTITE's protocols based upon ISO-10993 biocompatibility standards, *with certificates available on request*. These medical device adhesives are ideally suited for the assembly of a variety of devices including but not limited to:

- Catheters
- Blood collection sets
- Syringes & lancets
- Drug delivery devices
- Anesthesia, epidural & insulin syringes
- Blood oxygenators
- Blood filters
- Heat exchangers
- Blood bowls
- Manifolds
- Blood pressure transducers

KEY PROPERTIES*	AA 3961™	AA 3963™
Viscosity (cP)	60	300
Appearance (uncured)	Clear, colorless	Clear, colorless
Fluorescence (uncured/cured)	High	High
Fixture Speed with 405 NM LED Curing System (secs)	< 5	< 5
Shear Strength on Polycarbonate Blocks (PSI)	3,000	3,400
Needle Pullout Strength on Polycarbonate Hubs & 27 Gauge SS Needle (lb-f)	36	33
Shore Hardness (Shore D)	71	72

*Estimated values. Please consult the LOCTITE® TDS for additional data and details.

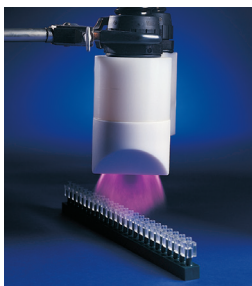

SURFACE TREATMENT:


Polyolefins such as polyethylene (PE) and polypropylene (PP) are often used for disposable medical devices due to their low cost and versatile properties. Since PE and PP are difficult-to-bond plastics, it is common to treat the plastic before assembly to increase the pull strength. Corona, plasma and other pretreatment methods have been shown to increase bond strengths on difficult-to-bond materials when using LOCTITE® AA 3961™ or AA 3963™.

ORDERING INFORMATION	AA 3961™	AA 3963™
25 ml Syringe	2464890	2483476
1 Liter Bottle	2464891	2483477



USA
 Henkel Corporation
 Engineering Adhesives
 One Henkel Way
 Rocky Hill, Connecticut 06067
 Tel: 1.800.LOCTITE (562.8483)
 Tel: 860.571.5100
 Fax: 860.571.5465
na.henkel-adhesives.com/medicaladhesives

CANADA
 Henkel Canada Corporation
 2515 Meadowpine Blvd.
 Mississauga, Ontario, Canada L5N 6C3
 Tel: 1.800.263.5043 (within Canada)
 Tel: 905.814.6511
 Fax: 905.814.5391

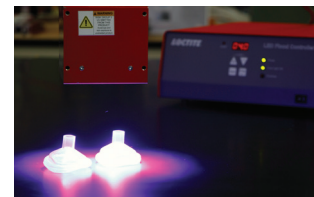


CURING

LOCTITE® AA 3961™ and AA 3963™ can be cured with UV or visible light. LOCTITE's newest LED curing units are ideally suited for rapid curing of select LOCTITE adhesives.



DISPENSING

LOCTITE® AA 3961™ and AA 3963™ can be dispensed manually or with LOCTITE® semi-automated and fully automated systems. See next.henkel-adhesives.com/us/en/products/all-shoppable-products.html for more details.

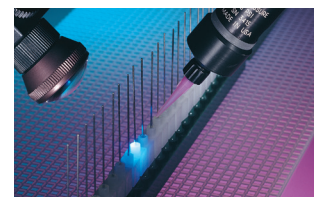

DETECTION

LOCTITE® AA 3961™ and AA 3963™ contain fluorescent additives allowing for detection in the uncured or cured state. Detection can be accomplished with a simple black light or with vision systems.

All marks used above are trademarks and/or registered trademarks of Henkel and its affiliates in the U.S. and elsewhere.
 ® = registered in the U.S. Patent and Trademark Office. © 2023 Henkel Corporation. All rights reserved. DSGN0030342 (06/23)

